

ADHESIVES

Choosing the wrong adhesive can leave you in a sticky situation. That's why ODP offers the popular BondPro™ adhesive system, complemented by a full line of bonding accessories. BondPro™ adhesives feature superior bonding characteristics ideal for metal, ceramic, and plastic/composite appliances. Because BondPro™ adhesives require no mixing, chair time is reduced, and the high tack paste minimizes bracket drift, allowing for quicker, easier and more accurate bracket placement. Stick with ODP's BondPro™ adhesive system, and you won't be left in a bind.

ONE STEP ORTHODONTIC ADHESIVE

Bond Pro™ One Step Orthodontic Adhesive

A unique composite resin system for bonding plastic, metal, or porcelain orthodontic brackets to etched enamel. The highly filled paste features a tacky composite viscosity that flows and eliminates bracket flotation. The Catalyst Resin/Paste system sets on contact, giving the clinician complete control. The primer also helps to seal the etched enamel. The cohesive strength of the adhesive is stronger than the bond to the tooth, allowing for one-piece removal after orthodontic treatment.



One Step Adhesive Kit

Kit Contents:

- 4x 5gm Syringes
- 12ml Primer
- 7ml Liquid Etchant
- Instructions
- 50 Sheet Mixing Pad
- 30x 2" Brushes
- 2 Mixing Spatulas

ORDER NUMBER M961K-01



Separate Refill Components

5gm Adhesive Syringe: M962-01
12ml Primer: M963-01
7ml Liquid Etchant: M964-01

One Step Adhesive Intro Kit

Kit Contents:

- 1x 5gm Syringe
- 3ml Liquid Etchant
- 5ml Primer
- Instructions

ORDER NUMBER M960K-01

LIGHT ACTIVATED ORTHODONTIC ADHESIVE

BondPro™ Light Activated Orthodontic Adhesive

A unique light-activated composite resin system for bonding plastic, metal, or porcelain orthodontic brackets to etched enamel. The paste features a tacky viscosity that flows well and eliminates bracket flotation. The light activated bonding resin/paste system sets on demand, giving the clinician complete control. The primer also helps to seal the etched enamel. The cohesive strength of the adhesive is stronger than the bond to the tooth, allowing for one-piece removal after orthodontic treatment with no enamel scarring. There is no mixing required.



Light Activated Adhesive Kit

Kit Contents:

- 4x 5gm Syringes
- 7ml Resin
- 7ml Liquid Etchant
- Instructions
- 50 Sheet Mixing Pad
- 30x 2" Brushes
- 2 Mixing Spatulas

ORDER NUMBER M971K-01



Light Activated Adhesive Intro Kit

Kit Contents:

- 1x 5gm Syringe
- 3ml Resin
- 3ml Liquid Etchant
- Instructions

ORDER NUMBER M970K-01

Separate Refill Components

5gm Adhesive Syringe: M972-01

7ml Resin: M973-01

7ml Liquid Etchant: M964-01

Adhesives

GLASS IONOMER BAND CEMENT

BondPro™ Automix Dual Cure Glass Ionomer Band Cement

A unique dual cure composite resin cement ideally suited for cementing metal molar bands. Chemically bonds to both tooth structure and metal band to create a strong anchor. The glass ionomer technology creates high tensile and compressive strength, reducing washout and loose bands. The cement is ideally mixed and has no strong odor or aftertaste. The adhesive can be light cured to immediately set band in place. Without light activation, a complete set is achieved in five minutes. Contains fluoride for added tooth protection. Easily dispensed from economical screw syringes or in smart automix system.



Automix Dual Cure Band Cement Kit

Kit Contents:

- 4x 10gm Syringes (5gm Catalyst & 5gm Base)
- 50 Sheet Mixing Pad
- 30 Dispensing Tips
- Instructions

ORDER NUMBER M976K-01

Separate Refill Components

10gm Adhesive Syringe: M977-01

GLASS IONOMER BAND CEMENT

Automix Dual Cure Band Cement Intro Kit

Kit Contents:

- 1x 10gm Syringe (5gm Catalyst & 5gm Base)
- 7 Dispensing Tips
- Instructions

ORDER NUMBER M975K-01



Glass Ionomer Powder/Liquid Band Cement

An advanced formulation designed as a multipurpose final cement for the cementation of orthodontic bands. The cement is radio-opaque and has cariostatic properties by naturally releasing fluoride.

Kit Contents:

- 18gm Powder
- 17.5ml Liquid
- 50 Sheet Mixing Pad
- Measuring Scoop
- Instructions

ORDER NUMBER M965K-01



Adhesives

BONDING SUPPLIES

BondPro™ Etchant Gel Kits

BondPro™ 37% phosphoric acid gel is a medium viscosity etchant gel used for etching enamel and dentin prior to bonding procedures. BondPro™ Etchant Gel is thixotropic, stays where placed and will not run. Does not dry out and has a long shelf life when properly stored. Washes off clean and is available in several syringe types.



Blue Etchant Gel Intro Kit

37% Phosphoric Acid

Kit Contents:

- 1x 12gm Syringe
- Disposable Mixing Tips
- Instructions

ORDER NUMBER M981K-01

BondPro™ Bond Enhancers

Applied prior to bonding for stronger bonds.

10cc Bottle.

Ordered Individually

Porcelain Enhancer: M955-01

Metal Enhancer: M956-01


Plastic Composite Enhancer: M957-01



Mini Curing Pen

Don't be in the dark when the power goes out! Our mini curing pen has a 100,000 hour solid state LED bulb with a lifetime warranty. Cures most composites in 20 to 40 seconds. Radiometer reading will not exceed 30. Cold sterilizable or autoclavable after batteries have been removed.

ORDER NUMBER M950-01

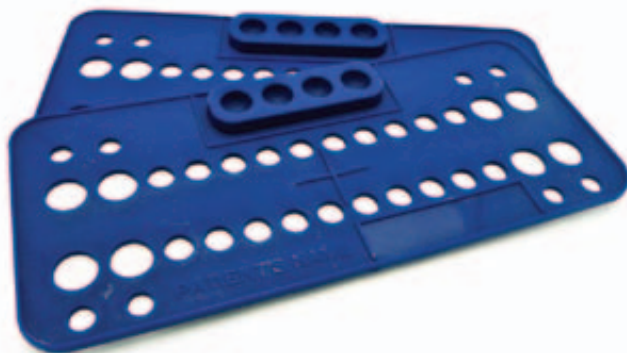

Plastic Bonding Trays

Plastic bonding trays are capable of holding one full bracket setup. Adhesive backing holds brackets securely in place. Available with mixing wells.

25Pk

ORDER NUMBER M990-02

BONDING RETRACTORS

Lip Retractors

Our transparent plastic lip retractors are ideal for use in taking photographs, or for difficult posterior bonding cases.

2Pk

MEDIUM: M984-01

SMALL: M985-01



Cheek Retractors

These flexible, comfortable retractors provide broad visibility and access for bonding. Autoclavable.

2Pk

MEDIUM: M986-01

SMALL: M987-01



Double Ended Photo Cheek Retractors

Durable plastic design is offered in 2 packs of autoclavable, heat sterilizable, or disposable. Enables maximum view of posterior while easily holding lips and cheeks.

2Pk

Autoclavable: M996-01

Heat Sterilizable: M997-01

Disposable: M998-01



Intraoral Cheek Retractors

The best retractors for bonding, offering unrestricted access to the entire mouth. Available in heat sterilizable, autoclavable, and 10/Pk disposable.

Ordered Individually

Autoclavable Small: M993-01

Autoclavable Medium: M993-02

Autoclavable Large: M993-03

10/Pk Disposable Small Green: M995-01

10/Pk Disposable Medium Blue: M995-02

10/Pk Disposable Large Red: M995-03